

UNITHERM WA 25

Orbital Shaking Waterbath



First Waterbath with powerful,
 noiseless magnetic drive.
 German technology

- Temperature range up to 80°C
- Speed range 10 - 220rpm
- Drive orbit 20mm
- Digital setting and programming time, temp. and speed and speed parameters
- Microprocessor control
- Overtemperature protection
- Acoustic alarm
- Stainless steel lid
- Different clamps available
- Special designed hoods available
- Stainless steel platform – removable

UNITHERM WA 25 with shaking platform and plastic clamps

Specifications:

Bath volume:	max. 25 liters								
Temp range:	5K over RT up to 80°C								
Shaking movement:	orbital								
Drive orbit:	20mm								
Speed range:	10...220rpm								
Heater consumption:	2000W								
Safety device:	overtemperature protection with acoustic alarm								
Dimensions (wxdxh):	<table> <tr> <td>exterior</td> <td>550x400x380mm</td> </tr> <tr> <td>interior</td> <td>380x330x230mm</td> </tr> <tr> <td>usable bath height</td> <td>160mm</td> </tr> <tr> <td>shaking platform</td> <td>350x300mm</td> </tr> </table>	exterior	550x400x380mm	interior	380x330x230mm	usable bath height	160mm	shaking platform	350x300mm
exterior	550x400x380mm								
interior	380x330x230mm								
usable bath height	160mm								
shaking platform	350x300mm								
Max. clamp capacity:	4x 1000ml, 6x 500ml, 8x 250ml, 16x 100ml, 30x 25ml;								
Weight, net:	approx. 28kg								
Power:	220V / 50Hz or 110V / 60Hz								


Dimensions and specifications are subject to change without notice for further improvements on design.
 Dimensions may vary dependent on different accessories be used.

CLAMPS IN DIFFERENT SIZES:

made of plastic:
100ml, 250ml, 500ml, 1000ml;



made of stainless steel:
10ml, 25ml, 50ml, 100ml, 250ml, 500ml, 1000ml,
2000ml, 3000ml, 5000ml;

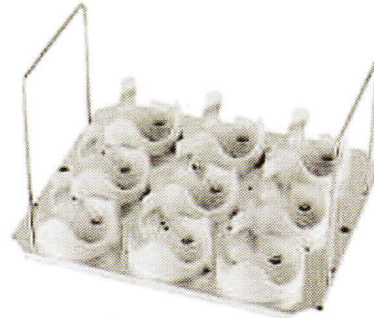

Standard platform with clamps for 250ml flasks

UNITHERM HYBRIDIZATION WATER BATH



- Temperature range up to 100°C
- Microprocessor controlled
- Light weight, fast heating water bath
- Combined with UNITWIST 3D shaker
- Speed range 10-100rpm
- UNITWIST spend power for water bath

UNITHERM HB Hybridization water bath + UNITWIST 3D

Dimensions and specifications are subject to change without notice for further improvements on design.
Dimensions may vary dependent on different accessories be used.